


# Fitting Guide

## Nathos Nova PR Slimtube

### Nathos™ Nova PR Slimtube

- 1 Right or left indicator\* on back of hearing aid
- 2 Charging contacts
- 3 Multi-function button with built-in indicator light
- 4 Slim tube
- 5 Retention cord
- 6 Speaker (with Dome earpiece attached)

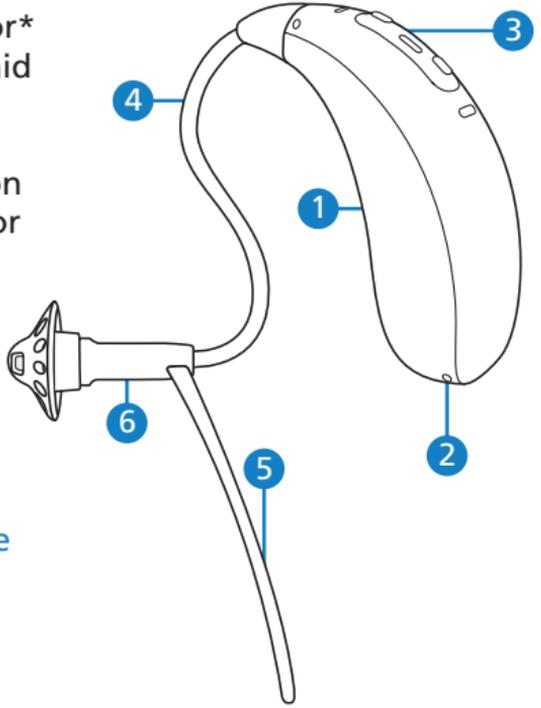

\*Right = red | Left = blue


### Turning on your hearing aid

Switch on/off the hearing system by pressing firmly on the lower part of the button for 3 seconds until the indicator light blinks

On: indicator light is solid green


Off: indicator light is solid red

*Please note:* Hearing aid automatically turns on when removed from the charger

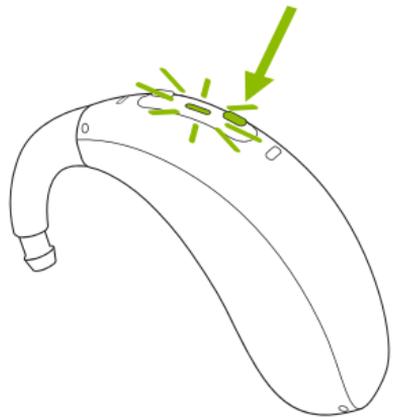

## Charging your hearing aid

Attach the UK adapter to the power lead, and plug this into a power outlet. Plug the smaller end of the cable into the USB-C port on the back of the charger. A indicator light on the charger, next to the cable will light up green when the charger is connected to the power supply.

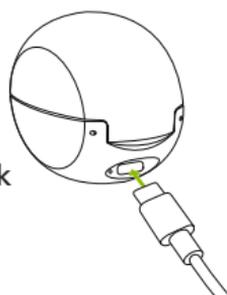

Insert the hearing aid into the charging ports. Make sure the left and right hearing aid markings match the left (blue), right (red) indicators on the charging ports.



The indicator light will show the charging state of the battery until the hearing aid is fully charged. When fully charged the indicator light will remain on, i.e. solid green light.



Remove the hearing aids from the charging inserts by

1. Pulling the hearing aids gently toward yourself and
2. Lifting them out of the charger.


The hearing aid switches on automatically when removed from the charger.



## Inserting the hearing aid



1 Place the hearing aid behind your ear.

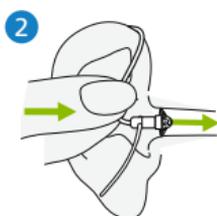

2 Insert the earpiece into your ear canal.



3 If there is a retention cord attached to the earpiece, tuck it into the bowl of your ear to secure your hearing aid.

## Changing the volume/program

The Multi-function button functions as a volume control and/or a program change depending on the hearing aid programming.

