


How to clean your open fit hearing aid and slim tube



1. Unscrew the slim tube to separate from the hearing aid.
2. Thread the cleaning wire through from the top to the bottom of the slim tube.
3. Pull the cleaning wire through the slim tube then clean wire with a damp cloth.
4. To re-attach the slim tube, push on to the hearing aid.

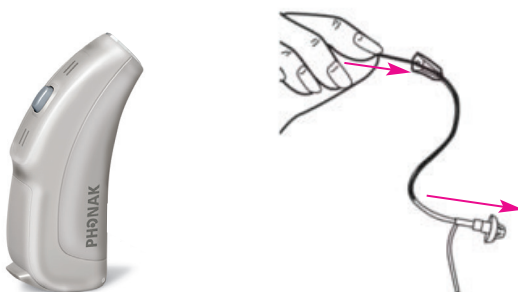


Audiology Patient Information

Quick Fit Guide - Slim Tube Fit


How to fit your open fit hearing aid and slim tube




Place the hearing aid over the top of your ear so that it sits behind your ear. The tip with dome needs to point towards your ear canal.

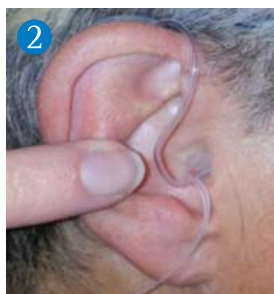

Whilst leaving the hearing aid behind your ear, push the tip with dome into your ear canal.

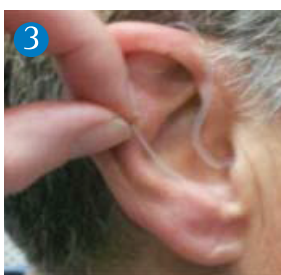

Once this is in, place the retaining clip into the bowl of your ear.

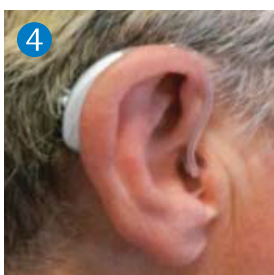

The tip with dome sits in the ear canal, the retainer in the bowl of the ear and the hearing aid is tucked behind the ear, correctly shown in this picture.